




Swiftstyle
specialist equipment




A fresh approach to Postural Stability!

We have put a lot of effort into designing products that work well – especially with our large client group that have very strong tone patterns. We believe that the Postureflex range rises to these challenges.

EASY TO USE – whether you are fitting the harnesses to a wheelchair, where the straps go over the backrest, or a buggy, where they fit through the backrest, you will find the securing the harness makes sense every time. That's because we have been in your situation too.

A BIT OF MOVEMENT – we probably should call our harnesses 'dynamic'. What we mean is that they have a little bit of give. This characteristic is normally achieved using Neoprene, - but we use a Heavy Duty Spacer Fabric – in two layers as you see below. This is an engineered fabric, which gives a realistic limit to the elasticity.

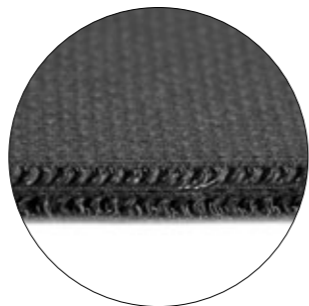

BREATHABILITY – one of the many benefits of steering away from Neoprene - our Spacer Fabric allows the heat and moisture to escape. This is critical for long term wheelchair users

POSTUREFLEX[®]



Using Harnesses to achieve Postural Alignment

1st Golden Rule - Start with the Hips. We want to align the pelvis which depends on the Ischial Tuberosities being located comfortably at the rear of the seat cushion. An Ischial ramp on the seat will help to achieve an anterior tilt of the pelvis

To maintain this pelvis position, the Hip belt should come as low as possible, locating on the pelvic bone. A 4 point hip belt has front pull anchor straps that prevent the belt from riding up into the abdominal cavity. So now we have the Pelvis central and parallel with the back of the seat.

UPPER BODY HARNESS

We aim to minimise the level of support when stabilising the upper body in order to make the user as comfortable as possible. We have a range of options that offer support in different areas.



SECURING FEET

We often find that if feet are not secured, users with a lot of leg movement very quickly lose their correct position at the hips. In addition to this, there are those hazards of feet slipping off the side or the back of the footplates



Hip Belts

CHOOSING THE STYLE – We use two different types of Buckle – the Plastic Side Release Buckle that gives a dual pull tensioning adjustment and the Centre Release Metal buckle which has one tensioning strap. You can also either choose a 2 point belt as shown in the top picture or a 4 point hip belt as shown in the lower two pictures.

CHOOSING THE SIZE – refer to this table to select the appropriate size hip belt

Size	Hip Width (mm)	Hip Width (ins)
X-Small	200-300	8-12
Small	275-375	11-15
Medium	350-450	14-18
Large	425-550	17-22

PRODUCT SELECTION - 2 POINT HIP BELT

Size	Side Release (Plastic Buckle)	Price (Each)	Centre Release (Metal Buckle)	Price (Each)
X-Small	HBX2D	£49.50	HBX2C	£64.50
Small	HBS2D	£49.50	HBS2C	£64.50
Medium	HBM2D	£53.50	HBM2C	£69.50
Large	HBL2D	£53.50	HBL2C	£69.50

PRODUCT SELECTION - 4 POINT HIP BELT

Size	Side Release (Plastic Buckle)	Price (Each)	Centre Release (Metal Buckle)	Price (Each)
X-Small	HBX4D	£59.50	HBX4C	£69.50
Small	HBS4D	£59.50	HBS4C	£69.50
Medium	HBM4D	£62.50	HBM4C	£72.50
Large	HBL4D	£62.50	HBL4C	£72.50



Foot Straps

CHOOSING THE STYLE – We have two different types of Foot Fixing. Our top picture shows the 4 point Foot Fast, which is designed to work with standard Remploy footplates. This is secured with a single Velcro strap that pulls the foot backwards and downwards at the same time. This one is particularly useful for users that have rotation in the lower leg, as it more successfully maintains foot alignment. The 4 point Foot Fast also remains open while not in use, making it very easy to secure the foot.

Our Lower Picture shows our conventional Footfasts that hug the ankle. These are very comfortable, and can be used in conjunction with toe straps if required.

FOOT FASTS

Size	Code	Price/Pair
Small	FF801	£79.50
Medium	FF802	£79.50
Large	FF803	£79.50

4 POINT FOOT FASTS

Size	Code	Price/Pair
Small	FF4901	£84.50
Medium	FF4902	£84.50
Large	FF4903	£84.50





Upper Body Harnessing

Size	Chest Width	Shoulder Height
X-Small	125-200mm	280-350mm
Small	200-300mm	330-430mm
Medium	280-350mm	400-580mm
Large	330-450mm	450-630mm



The Centreline harness is a very comfortable harness designed to follow the upper body contours, with a good position for feeding tubes if required.

Size	Code	Price/Each
Small	CHH401	£89.50
Medium	CHH402	£89.50
Large	CHH403	£95.50



The conventional H-Style Harness provides a minimalist option with very good adjustment

Size	Code	Price/Each
Small	CHH501	£74.50
Medium	CHH502	£74.50
Large	CHH503	£79.50



SUPPORT WHERE IT'S NEEDED

The vest Harness gives more support across the chest where this is needed.

Size	Code	Price/Each
Small	VH601	£84.50
Medium	VH602	£84.50
Large	VH603	£89.50



The Six Point Vest Harness gives extra lateral support, with tension adjustability helping you to correct specific postural needs

Size	Code	Price/Each
X-Small	FF700	£85.50
Small	FF701	£89.50
Medium	FF702	£89.50
Large	FF703	£95.50

Accessories

All harnesses are supplied with Tube Clamp Fixing kits to enable you to attach the harnesses to various pieces of equipment



Swiftstyle Supplies

Delta House, Belper Road,
Heaton Mersey Industrial Estate,
Stockport SK4 3QW

Tel: 0161 432 0066

Fax: 0844 8000 232

Email: sales@swiftstylesupplies.co.uk


Swiftstyle
specialist equipment